





Solutions for Engineering & Manufacturing



HTI Cupping prototype UI

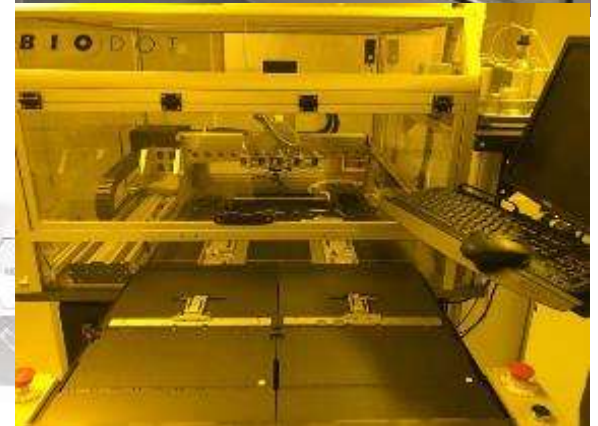


HEALTHCARE TECHNOLOGY INTL. LTD.

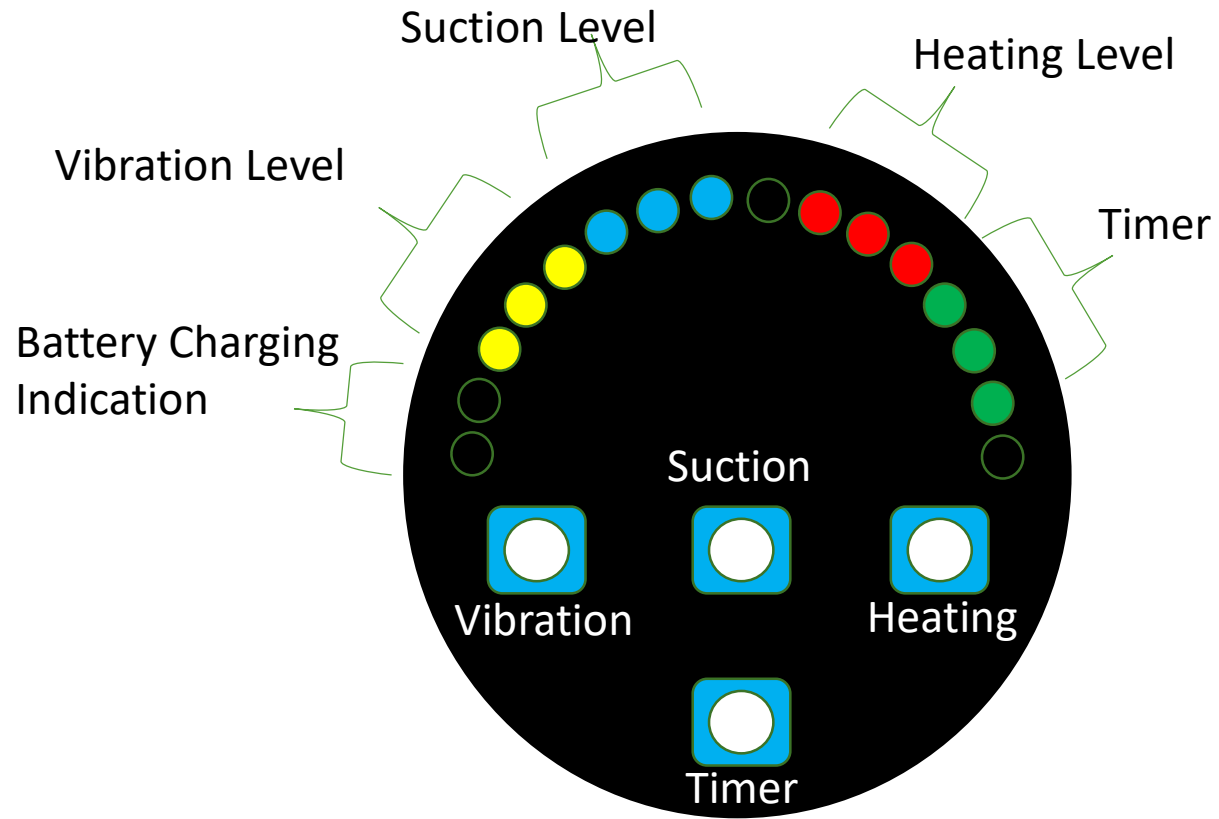

Vibration Level

- Press “Vibration” button
- 3 levels with PWM controlled vibration (100%,60%,30%)
- number of LED lighted represented the level strength


Vibration Level

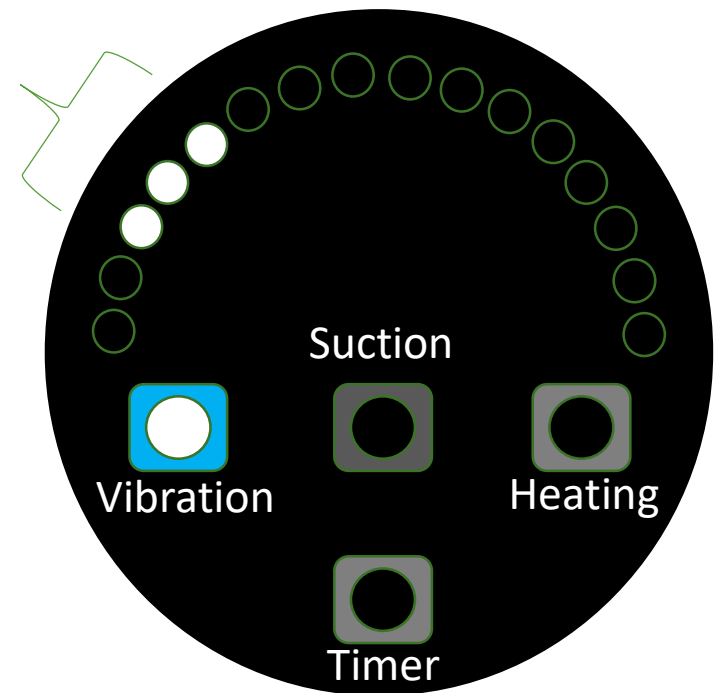

Suction Level


- Press “Suction” button
- 3 levels of suction pressure maintained by means of PWM control
- number of LED lighted defines the negative pressure (TBD)

Suction Level

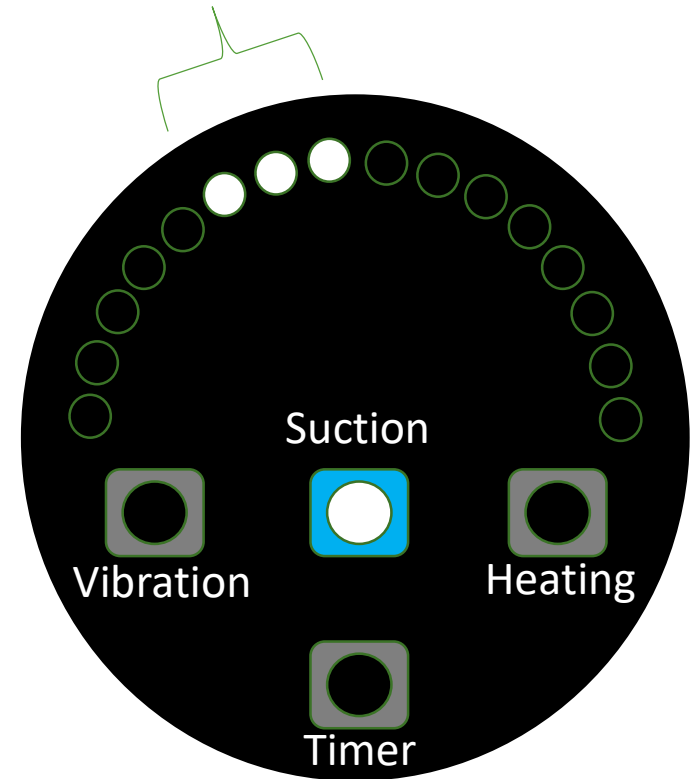

Heating Level

- Press “Heating” button
- 3 temperature of maintained by means of PWM control
- number of LED lighted represented the 3 selected temperature (assume 37,40,45)

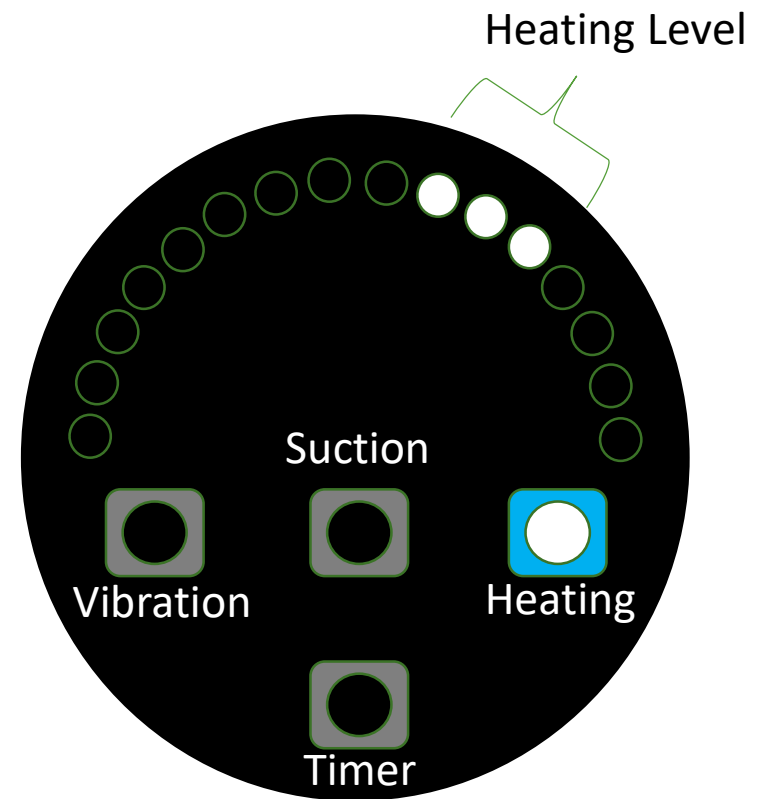

Timer Level

- Press “Timer” button
- 3 treatment time available
- number of LED lighted represent 5min,10 min or 15 min treatment time

